

# CP-500

## Chart Projector



## About Us

### Rex + Max = Rexxam

Rexxam, which means 'the king of the kings', is a respected and reliable brand.

Rexxam is a Japanese company with a celebrated 60 year history. With over 3,000 employees worldwide, Rexxam manufacture a wide range of products for various industries; from factory automation, automobiles and air conditioning systems, to beer and ski boots.

Since 1986, Rexxam has manufactured various high quality products for leading brands in the eye care industry, including SHIN-NIPPON. Rexxam had developed and manufactured products for SHIN-NIPPON since 1993 and in 2014 the company took over the SHIN-NIPPON brand.

We will be bringing high quality ophthalmic equipment to a global market. By combining precision engineering with industry leading innovation and experience in mass production, Rexxam produce unique products to support eye care specialists across the world.

Quality in vision care, we are Rexxam.



1960  
Foundation of Rexxam

1986  
Rexxam started the development and manufacturing of ophthalmic devices as an OEM supplier

1993  
Rexxam became the main OEM partner for SHIN-NIPPON  
**SHIN-NIPPON**

2014  
Rexxam acquired the SHIN-NIPPON brand  
**SHIN-NIPPON** by **Rexxam**

2018  
The manufacturer brand Rexxam was inaugurated  
**Rexxam**

**Rexxam**  
Quality in vision care

Proudly  Made in Japan



# LED Chart Projector

Rexxam CP-500 is a highly versatile LED chart projector that ensures bright and homogeneous test chart projection. Low energy consumption and no heat generation are guaranteed when using LED. LEDs have a long lifespan, which ensures years of uninterrupted use without having to replace the light source. You can select from A-Chart (Decimal Type) or B-Chart (Snellen ft Type). Projection distance ranges between 2m to 7m.

- White LED, adjustable to room environment
- Projection distance 2 to 7 metres
- Two models (Decimal Type and Snellen ft Type)
- Chart selection speed 0.3 seconds
- Fully integrated in Rexxam Digital Refractor DR-900

Compact & Solid








## CP-500 Remote Controller



## 50 charts with various masks

You can choose from Type-A Chart (Decimal) or Type-B Chart (Snellen ft). Each type offers 50 charts with various masks. Essential charts such as Landolt, Snellen, Letter and Children are available with horizontal and vertical line masks, as well as a single letter mask.

Additionally, various Balance Test Charts and Binocular Vision Test Charts such as Stereo, Coincidence (H/V) and Schober.

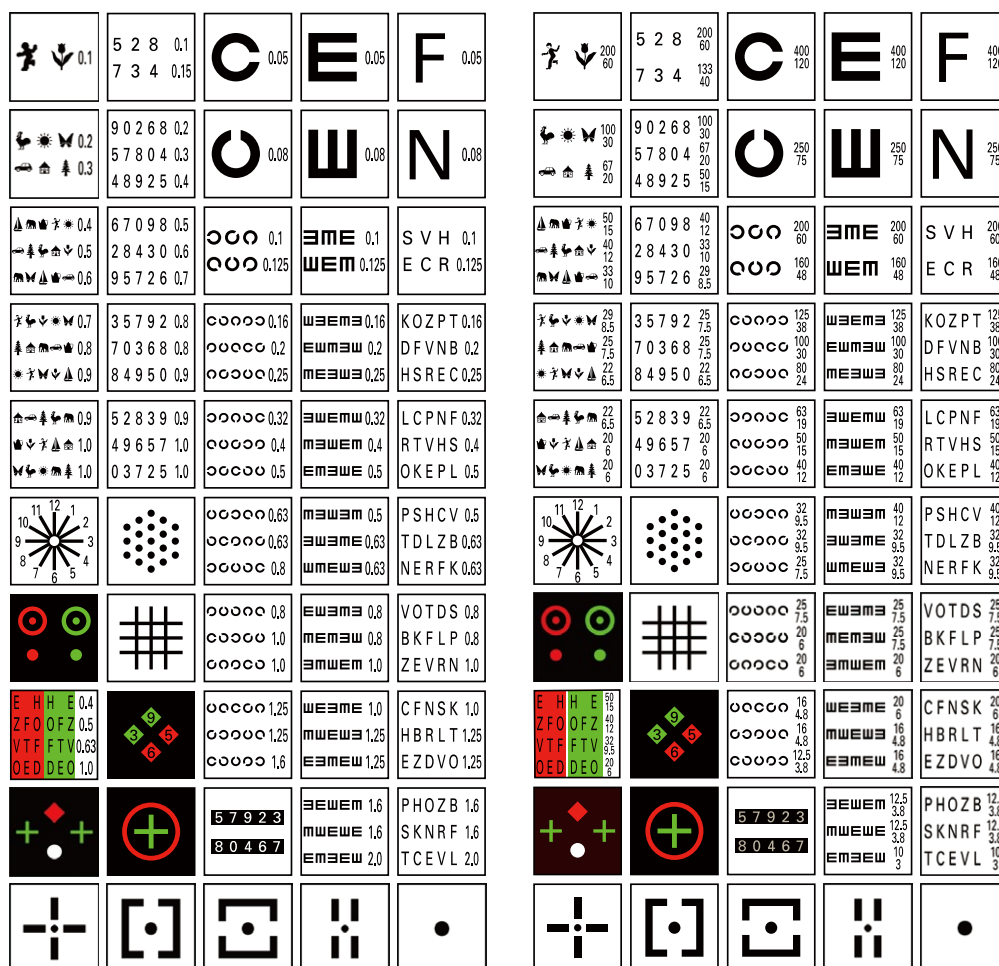

Type A Chart (Decimal)

Type B Chart (Snellen ft)



|                                |                   |  |
|--------------------------------|-------------------|--|
| Number of charts and masks     |                   | 50 charts and 38 masks                 |
| Luminance                      |                   | min.80 - max.320 Cd/m2(200Cd/m2 at 5m) |
| Chart and mask selection speed |                   | 0.3 seconds                            |
| Linear polarization            |                   | 45° ~ 135° , 135° ~ 45°                |
| Program                        |                   | 2 types (48 charts or masks)           |
| Tilting angle                  |                   | ±10°                                   |
| Auto power off                 |                   | After 5 min.                           |
| Lamp                           |                   | White LED                              |
| Projection distance            |                   | 2 ~ 7 m                                |
| Projection magnification       |                   | 30x (at 5m)                            |
| Size                           | Weight            | approx. 4.5kg                          |
|                                | Dimensions        | 248mm(W) × 336mm(D) × 234mm(H)         |
| Power                          | Power Voltage     | AC 100 - 240 V                         |
|                                | Power Consumption | 30VA                                   |

## Included Items

- Remote Controller
- Projection Screen
- Test Scale Sheet
- Spare Fuse
- Dust Cover



Projection Screen  
530mm(W) × 340mm(D)

## Optional Items

- Table Base Stand
- Wall Mount
- Chart Projector Arm
- Polarization Filter Glass
- Red / Green Glass



Table Base Stand

Design and specifications are subject to change without prior notice.

### Manufacturer



TOWA MEDICAL INSTRUMENTS Co.,Ltd.

1698-2 Nakano,  
Nakano-shi, Nagano, Japanese  
383-0013 Japan

### Contact



**Rexxam Co.,Ltd.**  
**Eye-care Instruments Sales Dept. Tokyo Office**  
2-4-2 Kandatsukasa-machi, Chiyoda-ku, Tokyo,  
101-0048 Japan  
TEL:+81-3-6262-9471 FAX:+81-3-6262-9472  
E-mail:eye@rexexam.co.jp  
Website:https://www.rexexam.co.jp

Proudly Made in Japan

### Distributed by

